

Single cell transcriptomics


10x genomics Chromium

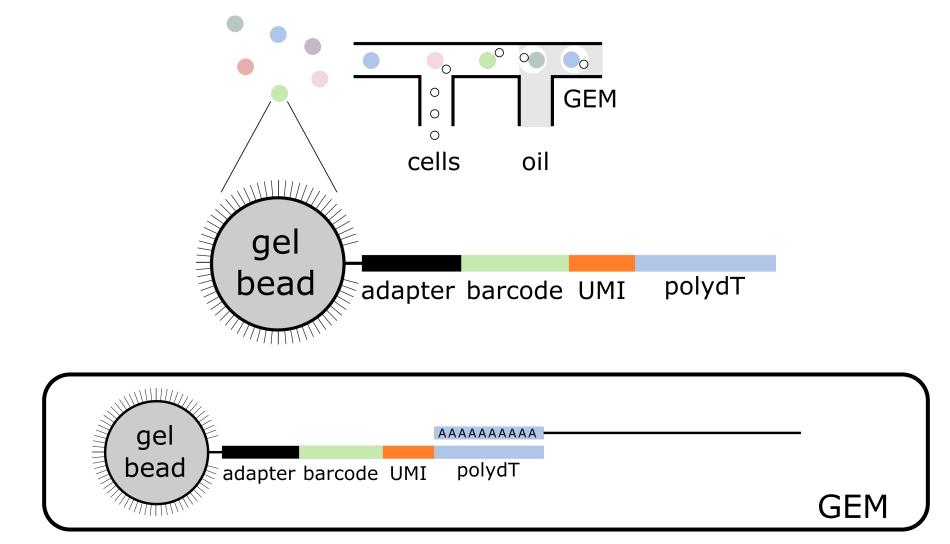

All captured **transcripts** from **single** cell: **identical** + **unique** barcode

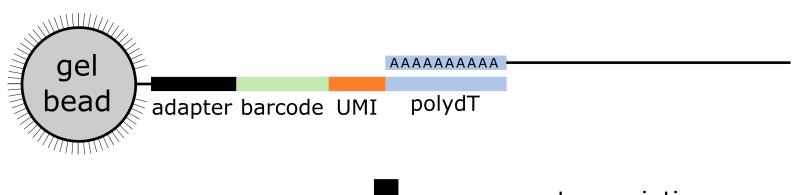

- reverse transcription
- breaking GEMs
- fragmentation
- primer ligation
- index PCR

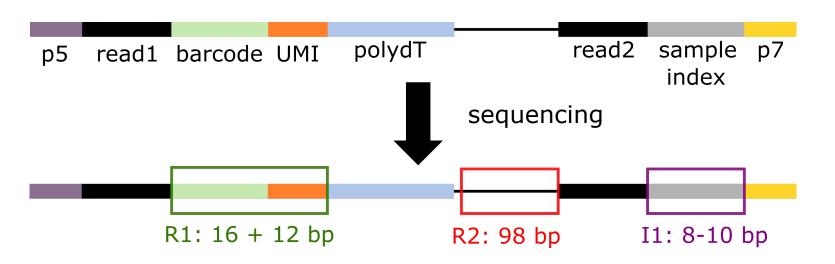

Sequencing output

```
ETV6-RUNX1_1_S1_L001_I1_001.fastq.gz
ETV6-RUNX1_1_S1_L001_R1_001.fastq.gz
ETV6-RUNX1_1_S1_L001_R2_001.fastq.gz
sample ID lane
```

- Dual indexing: second index in I2
- Indexes can also be added to fastq titles

After sequencing (preprocessing)

- 1. Assign reads to cell
- 2. Alignment
- 3. Quantification: # UMI/gene
- 4. Cell calling

For 10x all with cellranger count

Alternatives:

STARSolo Alevin

cellranger references

- Human & mouse: download pre-built from 10x website
- Other organisms: custom reference with cellranger mkref
- Exogenous marker genes (e.g. GFP): add sequence to both fasta and gtf
- Features (e.g.) hashing or surfaceproteins: feature barcode reference csv

Why count UMI (and not read alignments?)

- UMI: Unique Molecular Identifier:
 - Identifies each molecule (i.e. sequence) uniquely
- Molecules from a common PCR template
 - -> carry the same UMI
- By counting UMI: correct for PCR duplicates

Cellranger report

ETV6-RUNX1 1

Alerts

The analysis detected 🛕 1 warning.

	Alert	Value	Detail
A	Fraction of RNA read bases with Q-score >= 30 is low	59.4%	Fraction of RNA read bases with Q-score >= 30 should be above 65%. A lower fraction might indicate poor sequencing quality. This is Read 1 for the Single Cell 3' v1 chemistry and Single Cell 5' paired end, Read 2 for the Single Cell 3' v2/v3 chemistry and Single Cell 5' R2-only)
	30 13 tow		for the single cett 3 v2/v3 the mistry and single cett 3 v2-onty


Summary Analysis


3,091
Estimated Number of Cells

68,259 1,717

Mean Reads per Cell Median Genes per Cell


Sequencing ③			
Number of Reads	210,987,037		
Number of Short Reads Skipped	0		
Valid Barcodes	98.2%		
Valid UMIs	100.0%		
Sequencing Saturation	84.4%		
Q30 Bases in Barcode	96.4%		
Q30 Bases in RNA Read	59.4%		
Q30 Bases in UMI	96.5%		

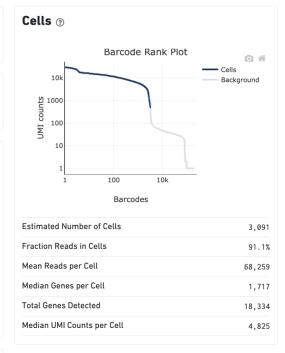


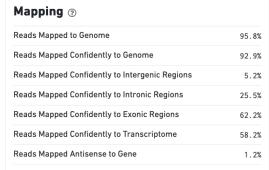

Sample				
Sample ID	ETV6-RUNX1_1			
Sample Description				
Chemistry	Single Cell 3' v2			
Include introns	False			
Reference Path	nger/refdata-cellranger-GRCh38-3.0.0			
Transcriptome	GRCh38-3.0.0			
Pipeline Version	cellranger-6.0.1			

Cell calling

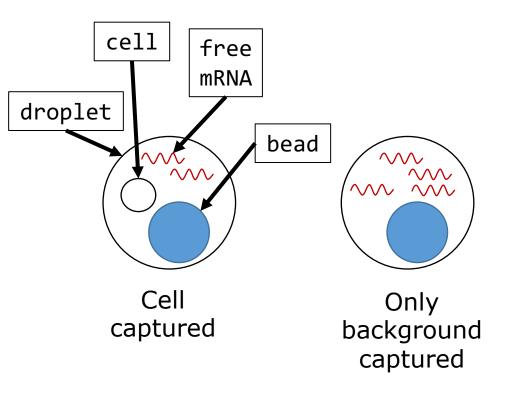

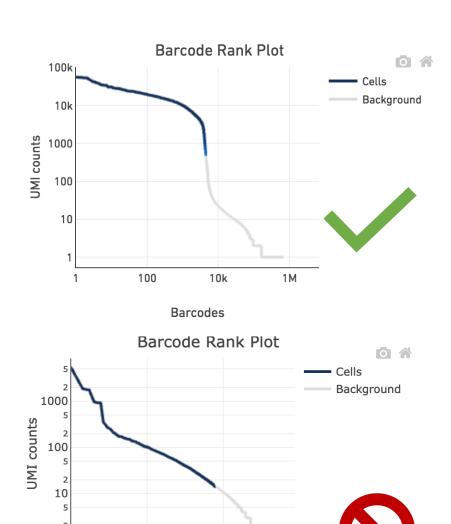

10k

100


Barcodes

Background 'cells': low #UMI/barcode

Other parameters

- Per channel:
 - Standard assay: 500-10,000
 - HT: 2,000-20,000 cells
 - LT: 100-1,000 cells
- Sequencing saturation
- Reads mapped to genome/transcriptome

$$saturation = 1 - \frac{\# unique \ reads}{\# \ reads}$$

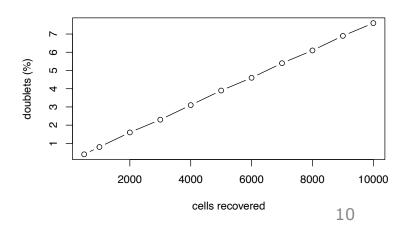

10x single cell flex

- FFPE fixed cells
- Based on probe hybridzation:
 - Specificy through ligation
 - ~3 probes/gene
 - Only human and mouse
 - Hybridized probes are sequenced
- 16 barcoded probe sets allows for multiplexing!

